

RESEARCH

Open Access



Post-thaw ATP supplementation enhances cryoprotective effect of iodixanol in rat spermatozoa

Suhee Kim, Sarah Hooper, Cansu Agca and Yuksel Agca*

Abstract

Background: Successful cryopreservation of rat spermatozoa from various strains still remains a challenge. The objective of this study was to determine if combinations of OptiPrep™ (iodixanol) and adenosine 5'-triphosphate (ATP) can improve rat sperm function during the cryopreservation procedure.

Methods: Epididymal rat spermatozoa were frozen under different OptiPrep™ concentrations (0, 1, 2, 3 or 4 %) and were diluted with media supplemented with or without 2 mM ATP after thawing. Post-thaw sperm motility, acrosomal membrane integrity (AMI) and mitochondrial membrane potential (MMP) were then evaluated. In addition, the effect of different OptiPrep™ concentrations on fresh and cooled rat spermatozoa was tested via motility.

Results: There was no effect of OptiPrep™ on motility of fresh and cooled spermatozoa. The supplementation of 1 and 2 % OptiPrep™ increased motility of frozen spermatozoa at 10 min after thawing, while it did not improve motility of spermatozoa at 3 h after thawing in the absence of ATP. During incubation of thawed spermatozoa, the ATP addition protected time-dependent decrease in motility after thawing in OptiPrep™-treated samples. OptiPrep™ had no effect on AMI and MMP in frozen-thawed spermatozoa but combinations of OptiPrep™ and ATP improved MMP in frozen-thawed spermatozoa.

Conclusions: Iodixanol has cryoprotective effects during rat sperm freezing without any toxic effect. Moreover, the combinations of iodixanol and ATP have a beneficial role in maintaining function of frozen-thawed rat spermatozoa for long period of incubation post-thaw.

Keywords: ATP, Cryopreservation, Rat spermatozoa, Iodixanol

Background

The rat has been recognized as one of the most valuable laboratory animals for biomedical and genomic research [1–4]. With the production of a large number of mutant or transgenic rat strains, there is an increasing demand for preservation and distribution of rat lines. Although rats are relatively easy to house and breed, it is expensive to maintain and transport genetically modified strains [5, 6]. For this reason, cryopreservation of gametes and embryos has been developed to maintain genetically valuable materials for an extended period and facilitate transportation of genetic materials at a low cost

[7, 8]. In particular, cryopreservation of spermatozoa has been well known as the most efficient way to preserve animal germplasm [9, 10], because of simplicity of epididymal sperm collection, convenience of cryopreservation method and economics for use of a few animals [11, 12]. However, intrauterine insemination (IUI) or in vitro fertilization (IVF) using cryopreserved rat spermatozoa are very limited by the low success of embryo developmental competence and pregnancy rate in most rat strains [10, 12]. Therefore, development of optimal protocols for rat sperm cryopreservation are currently ongoing [12–16].

Iodixanol was originally developed as an x-ray contrast agent [17]. Today, it is also widely applied as a medium for density-gradient centrifugation and cushioned centrifugation technique due to its inert, nontoxic, nonionic

* Correspondence: agcay@missouri.edu
Department of Veterinary Pathobiology, College of Veterinary Medicine,
University of Missouri-Columbia, 1600 East Rollins Street, Room W191,
Columbia, MO 65211, USA

characteristics [18–22]. Interestingly, addition of iodixanol to the freezing extender helped maintain membrane integrity and functionality during the cryopreservation process of bovine spermatozoa [23]. It was demonstrated that addition of minor amounts of iodixanol affected the glass transition temperature of the solution and altered the structure of the growing ice crystals. Currently, the use of iodixanol for rat sperm cryopreservation has not been reported. Here, we hypothesized that low level of iodixanol might have cryoprotective properties without exerting toxic effects during rat sperm cryopreservation.

Frozen-thawed viable spermatozoa often cannot maintain viability as long as fresh viable spermatozoa. Thus, the energy production and supplementation are important to support motility of frozen-thawed spermatozoa for successful cryopreservation [14, 24]. Spermatozoa may attain access to eggs by mobilizing metabolic energy production in the form of adenosine triphosphate (ATP) to drive motility [25]. The freezing extender containing exogenous ATP improved the cryosurvival of rat epididymal spermatozoa [14]. As modification of the above reports, we thought that ATP can extend lifespan of frozen spermatozoa post-thaw by supplying energy to thawed spermatozoa when frozen spermatozoa are incubated with thawing media containing ATP. Therefore, the aim of this study was to evaluate protective effect of iodixanol on rat spermatozoa during cryopreservation and to determine maintainable effect of ATP on function of frozen rat spermatozoa during incubation after thawing. The goal of this study was to develop an optimized cryopreservation condition for rat spermatozoa using iodixanol and ATP by offering proper condition for enhanced viability of rat spermatozoa after thawing as well as freezing.

Methods

Animals

Six sexually mature (12–15 weeks old) male rats (outbred Sprague-Dawley [SD] strain) were used as sperm donors. Rats were housed in conventional rat cages at 20–25 °C in a controlled light environment (10 h dark/14 h light) and provided free access to water and standard rodent chow. The rats were housed in accordance with the policies of the University of Missouri Animal Care and Use Committee and the Guide for the Care and Use of Laboratory Animals.

Experimental design

The effect of iodixanol (OptiPrep™; 60 % w/v iodixanol in water) was evaluated during cryopreservation of rat spermatozoa. The OptiPrep™ concentrations for this study were 0 % [control], 1, 2, 3 and 4 %. Epididymal spermatozoa was divided into three as following: 1) fresh spermatozoa were incubated with different concentrations

of OptiPrep™ for 1 h to test potential toxic effect of iodixanol, 2) spermatozoa were cooled in freezing media containing different concentrations of OptiPrep™ for 1 h and 3 h to determine iodixanol effect on cooling time, 3) spermatozoa were cooled and frozen in freezing media containing different concentrations of OptiPrep™ to evaluate iodixanol effect for cryopreservation.

To check the effect of ATP after thawing, frozen spermatozoa was incubated in media containing 0 mM or 2 mM ATP after thawing. Incubation time was 10 min and 3 h to evaluate the maintaining effect of ATP on sperm function for long incubation. Spermatozoa frozen with different concentrations of OptiPrep™ were thawed and incubated with ATP to evaluate synergistic effect of OptiPrep™ and ATP.

Sperm collection, freezing and thawing

All chemicals were purchased from Sigma Chemical (St. Louis, MO, USA) unless otherwise stated. Male rats were euthanized by CO₂ inhalation. The cauda epididymides were excised and placed in 35-mm petri dishes containing HEPES-buffered Tyrode lactate (TL-HEPES) media [26] supplemented with 3 mg/mL bovine serum albumin (BSA) and 0.11 mg/mL pyruvic acid. The cauda epididymides were dissected with a fine scissors to allow spermatozoa to swim out for 10 min at 37 °C. The sperm suspension (fresh spermatozoa) was gently drawn into a plastic Samco transfer pipette (San Fernando, CA, USA) with 2-mm inner diameter and placed in 1.5-mL Eppendorf tube. The initial sperm concentrations were approximately 50×10^6 spermatozoa/mL and only semen with ≥ 70 % motility was used for this study. The collected spermatozoa were processed individually.

The freezing extender used in this study was the TEST (TES-Tris buffer)-sucrose-egg yolk freezing extender comprised of 25.14 g/L N-tris (hydroxymethyl) methyl-2-aminoethanesulfonic acid (TES), 13.14 g/L Tris base, 27.34 g/L sucrose, 20 % (vol/vol) egg yolk and 0.75 % (vol/vol) Equex-Paste (Minitüb, Tiefenbach, Germany). The spermatozoa were frozen as previously described [27]. Briefly, one part of fresh spermatozoa was mixed with four parts of freezing extender. The diluted spermatozoa were cooled at 4 °C for 60 min and were loaded into 0.25-mL straws. The straws were placed in liquid nitrogen (LN₂) vapor for 10 min, plunged into LN₂ and stored at LN₂ for at least 7 days. The straws were thawed by holding them in water (37 °C) for 10 s. The 200 µL of thawed spermatozoa were then diluted with 800 µL of HEPES-buffered fertilization media (modified Krebs Ringer bicarbonate [mKRB]) at the ratio of 1:4 and incubated at 37 °C for evaluation.

Iodixanol treatment

Fresh spermatozoa were evaluated for motility after 1 h incubation at 37 °C in TL-HEPES media containing

different OptiPrep™ concentrations at the ratio of 1:4 or were kept at 4 °C after addition of freezing extender containing different OptiPrep™ concentration at the ratio of 1:4. The cooled samples were diluted with TL-HEPES (1:4 ratio) 1 and 3 h after cooling initiation, incubated at 37 °C for 10 min, and evaluated for motility. Meanwhile, spermatozoa cooled for 1 h were frozen, diluted with fertilization media (1:4 ratio) after thawing, and incubated for 10 min and 3 h at 37 °C for evaluation of motility. Thawed spermatozoa were also evaluated for acrosomal membrane integrity (AMI) and mitochondrial membrane potential (MMP) at 10 min after thawing.

ATP treatment

To determine if ATP can play a role as an energy source after thawing, frozen-thawed rat spermatozoa were diluted with fertilization media (1:4 ratio) containing ATP (0 or 2 mM). The samples were processed as described in above section (iodixanol treatment).

Evaluation of sperm function

Sperm motility

Twenty μL of sperm suspension were loaded into 80- μm deep dual-sided chamber (2 \times CELL; Hamilton Thorne Biosciences) at 37 °C and covered with cover glass. Computer-assisted sperm analysis (CASA, Hamilton Thorne Biosciences Inc, Beverly, MA, USA) system was then used to analyze the motility of rat spermatozoa placed onto the chamber [28]. Total and progressive motilities and average path velocity (VAP) were evaluated in six fields using CASA. The setting parameters and the definition of measured sperm motion parameters for the CASA were: frames per second, 20; duration of tracking time, 0.7 s; medium VAP cutoff, 25.0 $\mu\text{m}/\text{s}$; low VAP cutoff, 5.0 $\mu\text{m}/\text{s}$; count slow as motile, yes.

Sperm acrosomal membrane integrity

A PNA-Alexa Fluor 488/propidium iodide (PI) stain (L21409 and L7011, Molecular Probes Inc., Eugene, OR, USA, respectively) was used to determine AMI of rat spermatozoa. Aliquots of 200 μL of diluted frozen spermatozoa ($2 \sim 4 \times 10^6/\text{mL}$ spermatozoa) were mixed with 1 μL of 0.2 mg/mL PNA-Alexa Fluor 488 (final concentration, 1 $\mu\text{g}/\text{mL}$) and 1 μL of 200 μM PI (final concentration, 1 μM). The mixture was incubated at 37 °C for 15 min in 5 % CO_2 and then was analyzed by flow cytometry. Spermatozoa with PNA-/PI-, PNA+/PI- and PNA+ combined PNA+/PI- and PNA+/PI+ were considered to have viable spermatozoa with intact acrosomal membrane, viable spermatozoa with reacted acrosomal membrane and total spermatozoa with reacted/ or damaged acrosomal membrane, respectively.

Sperm mitochondrial membrane potential (MMP)

The sperm MMP was evaluated using the JC-1 fluorescent dye (Molecular Probes Inc.) as previously described [27]. Aliquots of 200 μL of diluted frozen spermatozoa ($2 \sim 4 \times 10^6$ spermatozoa/mL) were mixed with 1 μL of 100 μM JC-1 (final concentration, 0.5 μM) and 1 μL of 200 μM PI (final concentration, 2 μM) to exclude dead spermatozoa. The mixture was incubated at 37 °C for 30 min in 5 % CO_2 and then was analyzed by flow cytometry. The percentage of spermatozoa with viable and high MMP ($J_{\text{agg}}^+/\text{PI}^-$) and the Mean J_{agg} fluorescence intensity (J_{agg} MFI) per viable sperm were evaluated. To investigate the population of spermatozoa with high MMP of viable spermatozoa, the ratio of J_{agg}^+ to viable spermatozoa (%) was calculated as follows: $(J_{\text{agg}}^+/\text{PI}^- [\%] \div \text{PI}^- [\%]) \times 100$.

Flow cytometric analysis

All flow cytometry analyses were performed using a FACSCalibur flow cytometer (Becton Dickinson, San José, CA, USA) equipped with a 15 mW air-cooled 488 nm argon-ion laser and Cell Quest Pro software (Becton Dickinson). A total of 10,000 individual sperm-sized events were selected based on forward and side scatter and were collected at a flow rate of < 400 events/s. FL1 (PNA-Alexa Fluor 488), FL2 (J_{agg}) and FL3 (PI) signals were detected through 530 nm, 585 nm and > 670 nm band pass filter, respectively.

Statistical analysis

Statistical analysis was performed using SPSS software (version 17.0 for Windows; SPSS Inc., Chicago, IL). The Shapiro-Wilk test was utilized for normality analysis of the parameters. To compare the effect following OptiPrep™ concentrations, one-way repeated-measures (RM) analysis of variance (ANOVA) was used in data with a normal distribution and Bonferroni adjustment for multiple comparisons was used for *post hoc* analysis. Otherwise, the non-parametric Friedman test and *post-hoc* analysis with Wilcoxon signed ranks test was conducted with Bonferroni correction for data that was not normally distributed. For comparisons following incubation time and ATP supplementation, paired *t*-test or Wilcoxon signed ranks test was performed according to normality analysis. Statistical significance was set at $P < 0.05$ and all tests were two-tailed. Values were presented as the mean \pm standard error of the mean (SEM).

Results

The effect of iodixanol on fresh and cooled spermatozoa

When fresh spermatozoa were exposed to TL-HEPES with different concentrations of OptiPrep™, total motility, progressive motility and VAP were not different among samples ($P > 0.05$; Fig. 1). The motility of cooled

spermatozoa was also evaluated after 1 and 3 h cooling in different OptiPrep™ concentration (Fig. 2). There were no significant differences in motility of cooled spermatozoa among the different OptiPrep™ concentrations ($P > 0.05$). OptiPrep™ supplementation did not prevent total and progressive motility from decreasing according to cooling time (Fig. 2a). Total motility decreased at 3 h cooling as compared to 1 h cooling and progressive motility also decreased at 3 h cooling relative to 1 h cooling ($P < 0.05$) except for 2 % OptiPrep™ concentration. Duration of cooling had no significant effect on VAP (Fig. 2b). In conclusion, supplementation of OptiPrep™ did not protect decrease in sperm motility caused by longer periods of cooling except for progressive motility shown in 2 % OptiPrep™ concentration. These data suggest that iodixanol had no protective effect for cooling process and no adverse effects on motility of fresh and cooled rat spermatozoa as well.

The effect of iodixanol and ATP on motility of frozen-thawed spermatozoa

Samples frozen with different OptiPrep™ concentrations were evaluated for motility after thawing and incubating in the presence or absence of ATP (Fig. 3). When the samples were incubated for 10 min in the absence of ATP post-thaw, the 1 and 2 % OptiPrep™ concentrations were superior to the other groups in total sperm motility ($P < 0.05$; Fig. 3a). However, there was no difference in total motility among OptiPrep™ concentrations at 3 h incubation under ATP-free media. Total sperm motility of 0, 1, 2 and 4 % OptiPrep™ showed remarkable reduction for 3 h incubation as compared to 10 min incubation, while only the 3 % OptiPrep™ appeared to prevent reduction in total motility for long incubation time post-thaw in the absence of ATP (Fig. 3a).

Interestingly, total motility was maintained for 3 h incubation in presence of ATP, without any differences between 10 min and 3 h incubation in all OptiPrep™ concentrations (Fig. 3a). In particular, total motility of OptiPrep™-treated groups increased in ATP-supplemented

media as compared to ATP-free media at 3 h incubation (Fig. 3a). The incubation with ATP-supplemented media after thawing prevented decrease in rat sperm motility caused by longer incubation after freezing-thawing. The 1 % OptiPrep™ concentration showed a high percentage of total motility among different OptiPrep™ concentrations for 3 h incubation with ATP supplementation, indicating that combinations of iodixanol and ATP can be beneficial system to increase and maintain the number of motile rat spermatozoa even in longer incubation post-thaw.

There were no differences in progressive motility and VAP among OptiPrep™ concentrations at 10 min and 3 h incubation post-thaw (Fig. 3b and c). ATP improved VAP of 1 and 2 % OptiPrep™ concentrations at 10 min incubation post-thaw (Fig. 3c), supporting synergistic effect of iodixanol and ATP. Three h incubation decreased progressive sperm motility in control and 1 % OptiPrep™ and VAP in 1 % OptiPrep™ as compared to 10 min incubation post-thaw in absence of ATP, while it did not show a decrease in progressive sperm motility and VAP in all groups in the presence of ATP. In conclusion, low concentrations (1 and 2 %) of OptiPrep™ enhanced sperm motility with ATP-supplementation.

The effect of iodixanol and ATP on AMI of frozen-thawed spermatozoa

Sperm acrosomal status after freezing-thawing was divided into 4 categories depending on PNA/PI staining patterns using flow cytometry and epifluorescence microscope (Fig. 4a). There was no significant difference in the percentages of PNA-/PI- (viable and intact acrosome), PNA+/PI- (viable and reacted acrosome) and PNA+ (total reacted/or damaged acrosome) spermatozoa among samples frozen under different OptiPrep™ concentrations regardless of ATP presence (Fig. 4b-d). However, when the frozen-thawed samples were diluted with ATP-supplemented media, the percentage of PNA-/PI- spermatozoa increased in control, 1 and 2 % OptiPrep™ as compared to dilution with ATP-free

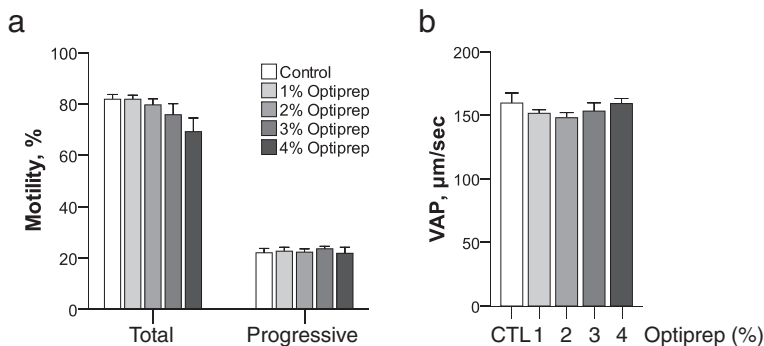


Fig. 1 The effect of OptiPrep™ on fresh rat sperm motility. The effect of OptiPrep™ on fresh rat sperm motility was evaluated using total and progressive motility (a) and VAP (b) after 1 h incubation with different concentrations of OptiPrep™ of fresh spermatozoa at 37 °C. $n = 6$

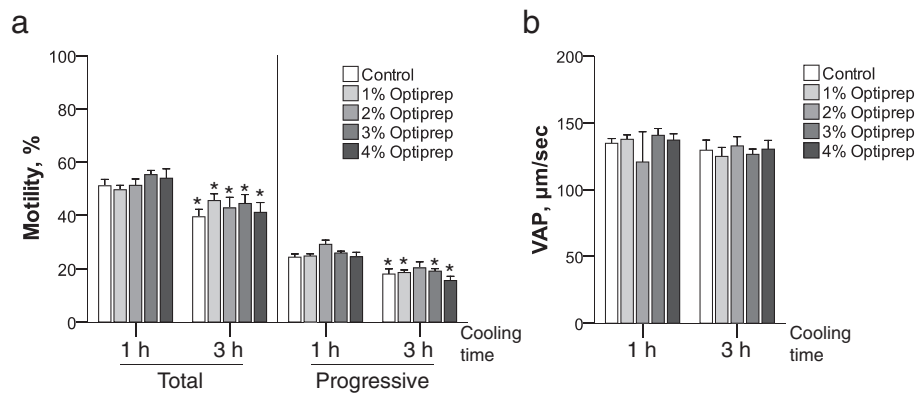


Fig. 2 The effect of OptiPrep™ on cooled rat sperm motility. The effect of OptiPrep™ on cooled rat sperm motility was evaluated using total and progressive motility (a) and VAP (b) 1 and 3 h after initiation of cooling to 4 °C with different concentrations of OptiPrep™ of rat spermatozoa. * $P < 0.05$ v.s. 1 h after cooling within same iodixanol treatment groups. $n = 6$

media (Fig. 4b). Meanwhile, the percentage of PNA+/PI− spermatozoa did not show differences in all groups by ATP treatment (Fig. 4c) and the percentage of PNA+ spermatozoa decreased in control, 1, 2 and 3 % OptiPrep™ by dilution with ATP-supplemented media as compared to dilution with ATP-free media (Fig. 4d). The OptiPrep™ had no effect on acrosome membrane while ATP can maintain stability of acrosome membrane in frozen-thawed rat spermatozoa.

The effect of iodixanol and ATP on MMP of frozen-thawed spermatozoa

Similar to effect on acrosome membrane of OptiPrep™ and ATP, there was no significant difference in MMP among samples frozen with different OptiPrep™ concentrations regardless of ATP (Fig. 4). However, ATP supplementation after thawing increased the percentage of Jagg+/PI− (viable and high MMP) spermatozoa and improved MMP of viable spermatozoa in groups treated with OptiPrep™ (Fig. 4e and f). Mean 75 % of viable frozen spermatozoa were Jagg + spermatozoa without significant differences following OptiPrep™ and ATP treatment (Fig. 4g). The combination of OptiPrep™ and ATP remarkably improved MMP function of frozen-thawed rat spermatozoa.


Discussion

Recently, it has been demonstrated that iodixanol is beneficial when added to a freezing extender for maintaining membrane integrity and functionality during the cryopreservation process of bovine and ram semen [23, 29]. Additionally, iodixanol contributed to preserving sperm motility after the freeze-thaw cycle. This study aimed to determine if iodixanol exhibited similar protective effects on rat spermatozoa during cryopreservation.

Supplementation of 5 % OptiPrep™ protected post-thaw progressive motility, membrane integrity, acrosomal

integrity and morphological damages in ram semen [29] and 2.5 % OptiPrep™ was superior for post-thaw motility in bovine spermatozoa [23]. In our study, low concentrations (1 or 2 %) of OptiPrep™ were superior to the higher concentration groups 10 min after thawing under ATP-free condition for total motility. The possible mechanisms that iodixanol protects spermatozoa during cryopreservation are attributed to an elevation of glass transition temperature and alteration of ice crystal formation in a non-colligative fashion, resulting in optimal environment for freezing [23]. Iodixanol alters ice crystal formation into intricate nets of dendrites in a concentrations-dependent manner during freezing [23]. Loose net of dendrites formed by lower concentrations of OptiPrep™ helped increase the cryoprotective effect on rat spermatozoa as compared to more compact net of dendrites formed by higher concentrations used in bovine and ram spermatozoa cryopreservation. Spermatozoa vary in size and shape according to species, which may be critical factors susceptible to structural changes of freezing condition by iodixanol. Falciform-shaped head and long tail of rat spermatozoa may serve as one of factors responsible for a better protective effect in lower concentrations of OptiPrep™ as compared to another species.

OptiPrep™ had a cryoprotective effect on freezing, not cooling, as shown in improved motility of 1 and 2 % OptiPrep™ at 10 min after thawing in our study. OptiPrep™ did not help preventing the decrease in total motility for longer storage of post-thaw as seen in 3 h incubation after thawing. Inevitable decrease in motility following time still occurred in cryopreserved rat spermatozoa in spite of OptiPrep™ treatment, with highly susceptible characteristics of rat spermatozoa to various external environments. In fact, when incubating for 3 h after thawing, sharp declines of total motility were observed even in 1 and 2 % OptiPrep™ groups, which showed an increase in total sperm motility at the early

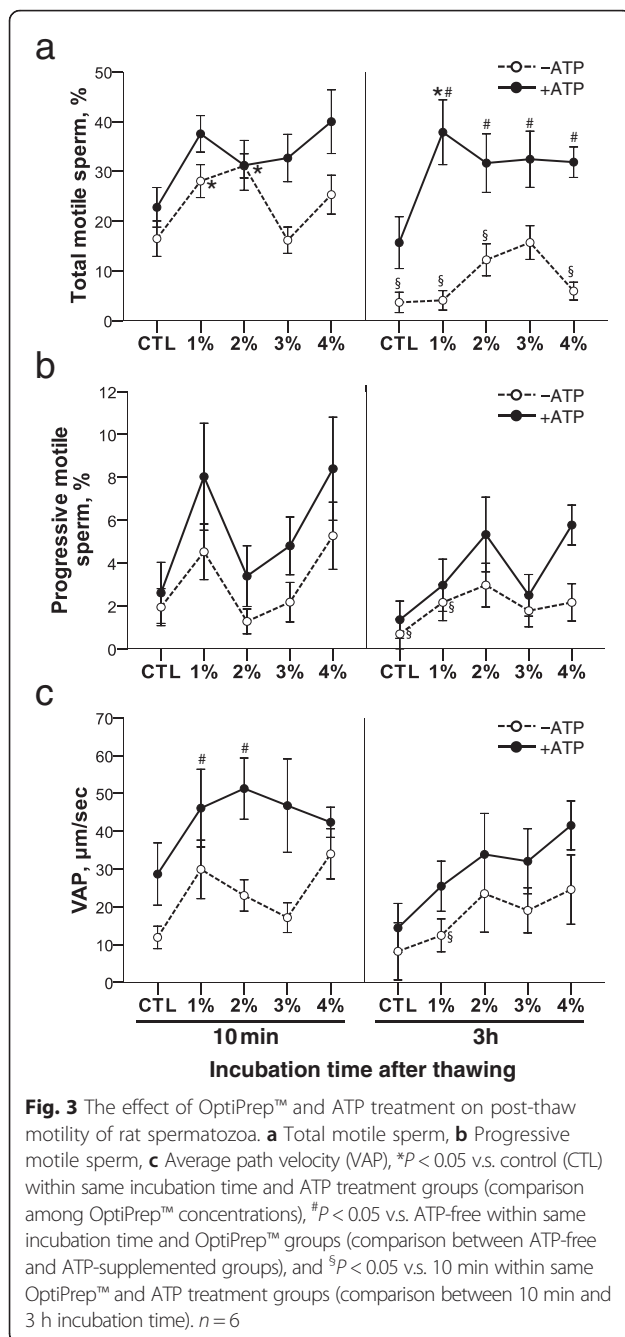

time point (10 min) after thawing in absence of ATP. Based on these results, long manipulation times consumed for AMI and MMP measurement may make difficult to determine the effect of OptiPrep™ on cryopreserved rat spermatozoa.

To overcome decreases in rat sperm motility caused by elapsed time after thawing, our second goal was to sustain motility of frozen-thawed rat spermatozoa by continuous supply of ATP. Both intracellular and extracellular ATP molecules play key roles in sperm function [30]. Intracellular ATP is the main energy source driving


sperm motility [30]. Extracellular ATP produces several downstream effects that improve sperm motility by increasing calcium levels or by activating cAMP signal transduction pathways [14]. However, there were no reports describing the effect of ATP on spermatozoa in thawing media as far as we know. This is first study to determine if supplement of ATP into thawing media can maintain the function of rat spermatozoa survived from cryodamage.

In this study, ATP supplementation to thawing medium prevented the decrease in total motility observed in groups incubated with ATP-free thawing medium for 3 h (Fig. 3a). ATP supplemented-OptiPrep™ groups showed improvement in total motility compared to ATP-free-OptiPrep™ groups at 3 h incubation of thawing (Fig. 3a). In particular, total motility of the 1 % OptiPrep™-treated group was superior to that of all other groups at 3 h incubation after thawing in ATP-containing media. This finding indicates that exogenous ATP can maintain sperm motility for longer periods of storage after thawing and combination of ATP and iodixanol can have a synergistic effect on enhancement of sperm motility. Therefore, ATP treatment with iodixanol may be beneficial during assisted fertilization by maintaining and extending activity of functional spermatozoa after freezing-thawing. In fact, extracellular ATP treatment of fresh mouse spermatozoa enhanced IVF success throughout alteration of motility by increased intracellular calcium level [30].

ATP supplementation also led to improvement in MMP in all groups treated with OptiPrep™ as compared to ATP-free condition. Moreover, the combination of ATP and OptiPrep™ had a strong synergistically effect in enhancing MMP function of frozen rat spermatozoa after thawing. However, considering that ATP increased AMI of frozen rat spermatozoa even in OptiPrep™-free group as well as 1 and 2 % OptiPrep™ groups, ATP seems to improve AMI of frozen rat spermatozoa alone regardless of OptiPrep™ treatment.

The effect of extracellular ATP treatment on sperm acrosome is controversial. Extracellular ATP was reported to increase acrosomal exocytosis (AE) in human and bovine fresh spermatozoa [31, 32] and did not affect AE in mouse fresh spermatozoa [30], while it increased the proportion of intact acrosomes during freezing as well as rat fresh sperm [14]. In the current study, ATP supplementation inhibited AE and increased acrosomal membrane integrity of cryopreserved rat spermatozoa in some OptiPrep™ groups as well as control group without OptiPrep™ treatment, indicating that ATP may affect capacity of rat spermatozoa to fertilize oocytes. The mechanisms that ATP influences on AMI will be studied in rat spermatozoa in the future.

Surgical intra uterine insemination (IUI) using rat sperm frozen in 2 % OptiPrep™ resulted in pregnancies and live fetuses (Additional file 1: Table S1). While 9 out

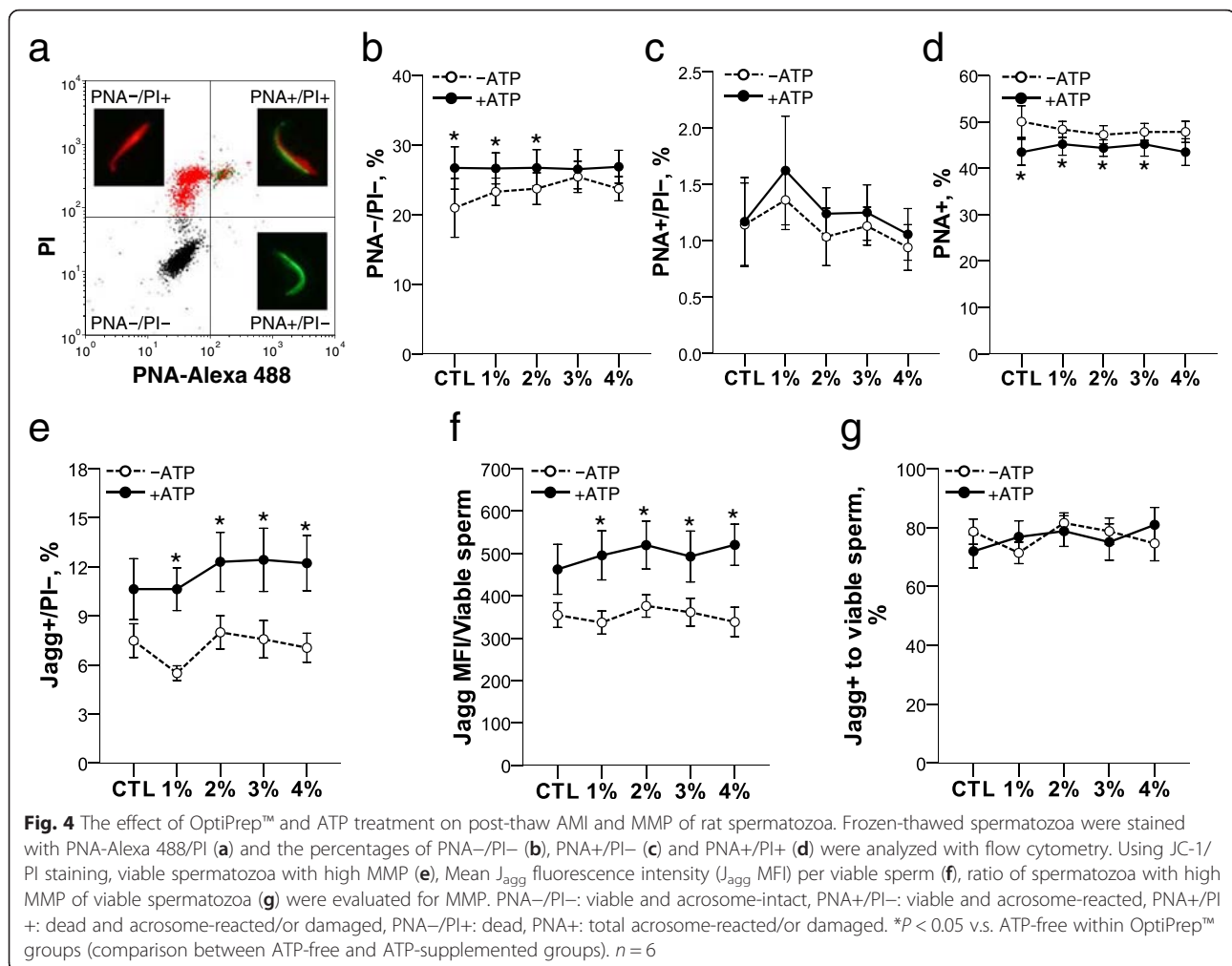

of 10 (90 %) rats become pregnant using fresh spermatozoa, 4 out of 11 (36 %) rats became pregnant using frozen-thaw sperm treated with 2 % OptiPrep™ (Additional file 1: Table S1). Although we did not compare the pregnancy and delivery rate between the absence and presence of OptiPrep™, this study led to a possibility for fertilization and pregnancy using frozen-thawed rat spermatozoa. Further testing on OptiPrep™ will be needed to overcome low fertilization and pregnancy rate of frozen rat spermatozoa [12, 33].

Conclusion

Iodixanol supplementation with low concentrations to freezing extender increased post-thaw motility. ATP protected spermatozoa from the decrease in motility by post-thaw time and improved AMI & MMP when combined with iodixanol. This study suggests that combinations of iodixanol and ATP may lead successful freezing and fertilization of rat spermatozoa by cryoprotective effect of iodixanol during freezing and maintenance ability of ATP after thawing.

Additional file

Additional file 1: The additional file included method and data for "Developmental competence of frozen-thaw rat spermatozoa in the presence of 2 % Optiprep™ after surgical intra uterine insemination." For more detail, see Additional file 1. (PDF 78 kb)

Abbreviations

AE: acrosomal exocytosis; AMI: acrosomal membrane integrity; ATP: adenosine 5'-triphosphate; BSA: bovine serum albumin; CASA: Computer-assisted sperm analysis; FL1: channel detecting green fluorescence; FL2: channel detecting orange fluorescence; FL3: channel detecting red fluorescence; IUI: intrauterine insemination; IVF: in vitro fertilization; Jagg: JC-1 aggregated form; Jagg+: highly MMP potential; JC-1: 5, 5', 6, 6'-tetrachloro-1, 1', 3, 3'-tetraethyl-benzimidazolylcarbocyanine iodide; MFI: mean fluorescence intensity; MMP: mitochondrial membrane potential; PI: propidium iodide; TES: N-tris (hydroxymethyl) methyl-2-aminoethanesulfonic acid.

Competing interests

The authors declare that they have no competing interests.

Authors' contributions

SK and YA conceived and designed the study. SK, SH and CA performed all experiments. SK analyzed, interpreted data and wrote the manuscript. All authors read and approved the final manuscript.

Acknowledgements

This work was supported by the NIH National Center for Research Resources grants R21 RR025913-02, P40 RR16939 and NIH T32OD011126.

Received: 23 November 2015 Accepted: 22 January 2016

References

- Jacob HJ, Kwitek AE. Rat genetics: attaching physiology and pharmacology to the genome. *Nat Rev Genet.* 2002;3(1):33–42.
- Zhou Q, Renard JP, Le Fric G, Brochard V, Beaujean N, Cherifi Y, et al. Generation of fertile cloned rats by regulating oocyte activation. *Science.* 2003;302(5648):1179.
- Gibbs RA, Weinstock GM, Metzker ML, Muzny DM, Sodergren EJ, Scherer S, et al. Genome sequence of the Brown Norway rat yields insights into mammalian evolution. *Nature.* 2004;428(6982):493–521.
- Lazar J, Moreno C, Jacob HJ, Kwitek AE. Impact of genomics on research in the rat. *Genome Res.* 2005;15(12):1717–28.
- Hagiwara M, Choi JH, Devireddy RV, Roberts KP, Wolkers WF, Makhlof A, et al. Cellular biophysics during freezing of rat and mouse sperm predicts post-thaw motility. *Biol Reprod.* 2009;81(4):700–6.
- Landel CP. Archiving mouse strains by cryopreservation. *Lab Anim (NY).* 2005;34(4):50–7.
- Woelders H, Matthijs A, Zuidberg CA, Chaveiro AE. Cryopreservation of boar semen: equilibrium freezing in the cryomicroscope and in straws. *Theriogenology.* 2005;63(2):383–95.
- Pinkert CA, Dyer TJ, Kooyman DL, Kiehm DJ. Characterization of transgenic livestock production. *Domest Anim Endocrinol.* 1990;7(1):1–18.
- Nakagata N. Use of cryopreservation techniques of embryos and spermatozoa for production of transgenic (Tg) mice and for maintenance of Tg mouse lines. *Lab Anim Sci.* 1996;46(2):236–8.
- Nakatsukasa E, Inomata T, Ikeda T, Shino M, Kashiwazaki N. Generation of live rat offspring by intrauterine insemination with epididymal spermatozoa cryopreserved at –196 degrees C. *Reproduction.* 2001;122(3):463–7.
- Nakagata N. Cryopreservation of mouse spermatozoa. *Mamm Genome.* 2000;11(7):572–6.
- Nakatsukasa E, Kashiwazaki N, Takizawa A, Shino M, Kitada K, Serikawa T, et al. Cryopreservation of spermatozoa from closed colonies, and inbred, spontaneous mutant, and transgenic strains of rats. *Comp Med.* 2003;53(6):639–41.
- Yamashiro H, Han YJ, Sugawara A, Tomioka I, Hoshino Y, Sato E. Freezability of rat epididymal sperm induced by raffinose in modified Krebs-Ringer bicarbonate (mKRB) based extender solution. *Cryobiology.* 2007;55(3):285–94.
- Yamashiro H, Toyomizu M, Kikusato M, Toyama N, Sugimura S, Hoshino Y, et al. Lactate and adenosine triphosphate in the extender enhance the cryosurvival of rat epididymal sperm. *J Am Assoc Lab Anim Sci.* 2010;49(2):160–6.
- Yamashiro H, Toyomizu M, Toyama N, Aono N, Sakurai M, Hiradate Y, et al. Extracellular ATP and dibutyl cAMP enhance the freezability of rat epididymal sperm. *J Am Assoc Lab Anim Sci.* 2010;49(2):167–72.
- Mazur P, Leibo SP, Seidel Jr GE. Cryopreservation of the germplasm of animals used in biological and medical research: importance, impact, status, and future directions. *Biol Reprod.* 2008;78(1):2–12.
- Kjeken R, Mousavi SA, Brech A, Gjoen T, Berg T. Fluid phase endocytosis of [¹²⁵I]iodixanol in rat liver parenchymal, endothelial and Kupffer cells. *Cell Tissue Res.* 2001;304(2):221–30.
- Flinck A, Gottfridsson B. Experiences with iohexol and iodixanol during cardioangiography in an unselected patient population. *Int J Cardiol.* 2001; 80(2–3):143–51.
- Matas C, Decuadro G, Martinez-Miro S, Gadea J. Evaluation of a cushioned method for centrifugation and processing for freezing boar semen. *Theriogenology.* 2007;67(5):1087–91.
- McCann CT, Chantler E. Properties of sperm separated using Percoll and IxaPrep density gradients. A comparison made using CASA, longevity, morphology and the acrosome reaction. *Int J Androl.* 2000;23(4):205–9.
- Saragusty J, Gacitua H, King R, Arav A. Post-mortem semen cryopreservation and characterization in two different endangered gazelle species (*Gazella gazella* and *Gazella dorcas*) and one subspecies (*Gazella gazelle acaiae*). *Theriogenology.* 2006;66(4):775–84.
- Saragusty J, Gacitua H, Pettit MT, Arav A. Directional freezing of equine semen in large volumes. *Reprod Domest Anim.* 2007;42(6):610–5.
- Saragusty J, Gacitua H, Rozenboim I, Arav A. Protective effects of iodixanol during bovine sperm cryopreservation. *Theriogenology.* 2009;71(9):1425–32.
- Long JA. Avian semen cryopreservation: what are the biological challenges? *Poult Sci.* 2006;85(2):232–6.
- Storey BT. Mammalian sperm metabolism: oxygen and sugar, friend and foe. *Int J Dev Biol.* 2008;52(5–6):427–37.
- Bavister BD, Leibfried ML, Lieberman G. Development of preimplantation embryos of the golden hamster in a defined culture medium. *Biol Reprod.* 1983;28(1):235–47.
- Kim S, Agca C, Agca Y. Changes in rat spermatozoa function after cooling, cryopreservation and centrifugation processes. *Cryobiology.* 2012;65(3):215–23.
- Varisli O, Uguz C, Agca C, Agca Y. Various physical stress factors on rat sperm motility, integrity of acrosome, and plasma membrane. *J Androl.* 2009;30(1):75–86.
- Cirt U, Bagis H, Demir K, Agca C, Pabuccuoglu S, Varisli O, et al. Comparison of cryoprotective effects of iodixanol, trehalose and cysteamine on ram semen. *Anim Reprod Sci.* 2013;139(1–4):38–44.
- Rodriguez-Miranda E, Buffone MG, Edwards SE, Ord TS, Lin K, Sammel MD, et al. Extracellular adenosine 5'-triphosphate alters motility and improves the fertilizing capability of mouse sperm. *Biol Reprod.* 2008;79(1):164–71.
- Luria A, Rubinstein S, Lax Y, Breitbart H. Extracellular adenosine triphosphate stimulates acrosomal exocytosis in bovine spermatozoa via P2 purinoreceptor. *Biol Reprod.* 2002;66(2):429–37.
- Foresta C, Rossato M, Di Virgilio F. Extracellular ATP is a trigger for the acrosome reaction in human spermatozoa. *J Biol Chem.* 1992;267(27): 19443–7.
- Nakata M, Okuda Y, Yamashita Y, Nakauchi C, Ito J, Kashiwazaki N. Successful production of offspring using cryopreserved sperm via nonsurgical artificial insemination in rats. *J Reprod Dev.* 2012;58(4):501–4.

Submit your next manuscript to BioMed Central and we will help you at every step:

- We accept pre-submission inquiries
- Our selector tool helps you to find the most relevant journal
- We provide round the clock customer support
- Convenient online submission
- Thorough peer review
- Inclusion in PubMed and all major indexing services
- Maximum visibility for your research

Submit your manuscript at
www.biomedcentral.com/submit

